## IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be taken, including the following:

### READ ALL INSTRUCTIONS BEFORE USING

**DANGER** — To reduce the risk of electrocution:

- DO NOT handle the plug with wet hands.
- DO NOT place or store the charger base where it can fall or be pulled into a tub, sink, or toilet.
- **DO NOT reach for a charger base that has fallen into water; unplug immediately. Do not reach into the water under any circumstances. Have the unit checked by an Authorized Service Center after it has been in water.**
- The unit is designed for use in a 120V AC, 60 Hz outlet. Never force the plug into the outlet.
- DO NOT submerge the charger base in water or any other liquids.

**WARNING** — To reduce the risk of burns, electrocution, fire, or injury:

- Close supervision is necessary when this product is used by, on, or near children or disabled persons.
- Never use while bathing.
- Keep the supply cord away from heated surfaces.
- Never operate the unit if it has a damaged plug, if it is not working properly, if it has been dropped or damaged, or if it has been immersed in water. Refer to an Authorized Service Center for repair.
- DO NOT excessively bend or pull the cord. Unplug the cord by grasping the plug — NOT THE CORD.
- Use this product only for its intended use as described in this manual. Only use attachments recommended by the manufacturer.
- This product is designed for teeth cleaning only — not for any other purpose.
- To prevent damaged bristles from breaking off while the unit is in use, replace your brush head if bristles are worn or bent.
- **Always unplug this product immediately after using.**
- **Never use while sleeping or drowsy.**
- **Never drop or insert any object into any opening or hose.**
- **All servicing should be performed by an authorized service representative.**
- **This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.**
- **Children should be supervised to ensure that they do not play with the appliance.**
- **The supply cord cannot be replaced. If the cord is damaged, the appliance should be scrapped.**

### Having Problems? We Have Solutions!

**The bristles stop oscillating or slow down during use.**

- You may be applying too much pressure against the teeth. Unlike a manual brush, the OSCILL8™ Rechargeable Toothbrush cleans effectively without strong pressure.

- If your brush head stopped working soon after you began using it, try using a different toothpaste and rinse brush head thoroughly after each use. Regular gel toothpastes provide excellent cleaning and extend the life of brush heads.

- Try cleaning the brush head thoroughly with warm water or mouthwash.

- Take the brush head off the handle and then turn on the unit. If the shaft on the handle is not moving up and down rapidly, you probably need to replace your brush head.

**The instrument will not turn on.**

- If the shaft on the handle is not moving up and down rapidly, battery may need to be charged. **Check to see if the charger is plugged into a live outlet**, or if the plug is controlled by a wall switch. The blue LED light should illuminate. Charge the unit for 16 hours and try again.

- If charging is still unsuccessful, the battery may need to be replaced.

**My unit is losing power.**

- First, try a new brush head. Also, restore power to the battery by discharging fully and recharging for 16 hours.

**The power handle is warm to the touch.**

- This is normal and is caused by the induction system used to charge the instrument's batteries.

**WARNING:** For the purposes of recharging the battery, only use the charging base provided with this appliance.

- Do not use outdoors or operate where aerosol(spray) products are being used or where oxygen is being administered.
- **This appliance contains batteries that are non-replaceable.**

## SAVE THESE INSTRUCTIONS


**MEDICAL WARNING:** For your safety:

- Contact your dentist prior to use if you have had oral surgery in the previous 12 months.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if the bleeding continues after two weeks of use.
- Be sure that only the bristles touch your teeth and gums. Damage to teeth and gums may occur if any part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the OSCILL8™, and should always exercise particular care with any oral hygiene device.
- Avoid using a brush with crushed or bent bristles. Damaged bristles may break off while the unit is in use. In order to prevent this, replace your brush head every three months, or sooner, if the bristles become worn or bent.
- The OSCILL8™ is a personal care unit and is intended for household use only.

**WARNING:** Changes or modifications to this unit not expressly approved by Conair will void the warranty. Servicing should always be referred to an authorized service center.

### Battery Disposal

1. To remove the battery for disposal or recycling, use screwdriver (crosshead) to release the screw at the bottom side of power handle, grip the upper side of power handle and slide out the internal bracket.
2. Remove the tape on the first piece of battery. Firmly tap the bottom of the handle on a hard surface in order to slide out the first piece of battery. Tap again to slide out the second piece of battery.
3. Now that the battery is exposed, bring the battery to the nearest collection center for disposal or call the toll-free RBRC™ information line at 1-800-8-BATTERY. RBRC™ is a trademark of the Rechargeable Battery Recycling Corporation. DO NOT throw the batteries into the normal municipal waste system.

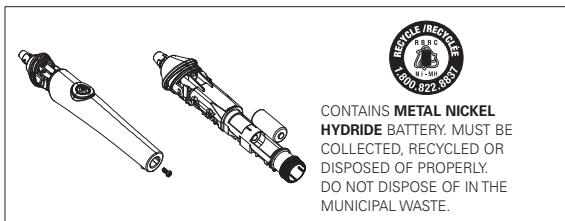CONTAINS METAL NICKEL HYDRIDE BATTERY. MUST BE COLLECTED, RECYCLED OR DISPOSED OF PROPERLY. DO NOT DISPOSE OF IN THE MUNICIPAL WASTE.

For more information on battery replacement for your power toothbrush, call 1-800-334-4031 (8am to 5pm Mountain Standard Time, Monday through Friday).

EVERY 6 MONTHS, COMPLETELY RUN DOWN AND RECHARGE YOUR POWER HANDLE TO MAXIMIZE THE LIFE OF THE RECHARGEABLE BATTERY.

Under various state and local laws, it is illegal to dispose of this battery in our municipal waste stream. The EPA-certified RBRC™ Battery Recycling Seal on the **nickel hydride (Ni-MH)** battery indicates Conair is voluntarily participating in an industry program to collect and recycle these batteries at the end of their useful life, when taken out of service in the United States or Canada. The RBRC™ program provides a convenient alternative to placing used **Ni-MH** batteries in the trash or the municipal waste stream, which may be illegal in your area. Please call 1-800-8-BATTERY for information on **Ni-MH** battery recycling and disposal bans/restrictions in your area. Conair's involvement in this program is part of our commitment to preserving our environment and conserving our natural resources.

RBRC™ is a trademark of the Rechargeable Battery Recycling Corporation.

### The OSCILL8™ Difference

The premium OSCILL8™ Rechargeable Toothbrush is designed to provide a high level of oral hygiene and improved oral healthcare at home. The revolutionary, 8-tuft brush head features independently oscillating tufts of tapered bristles that reach into areas around teeth and gums to eliminate plaque from even the tightest places. The compact brush head design comfortably fits into any size mouth, making this the ideal plaque remover for both children and adults.


### Using the OSCILL8™ Rechargeable Toothbrush: OSCILL8™ in 5 Easy Steps

#### 1. Charge the handle.

- Allow handle to charge for at least 16 hours before first use.
- Plug the charger base into an outlet, place the power handle on the base, and let it charge for 16 hours.
- The ring of blue light around the power switch will light up while charging.
- When fully charged, the ring light will go off. The LED light has an auto-off function after 16 hours.


#### 2. Attach the brush head.

- Attaches and removes with just a simple twist.

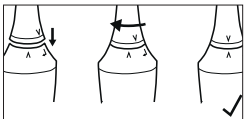- Insert the brush head with **V** mark pointing to **J** mark on power handle.
- Twist clockwise, until **V** mark on brush head and **A** mark on power handle are aligned.

#### 3. Wet the bristles, apply a pea-sized amount of toothpaste, and put in mouth.
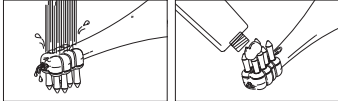
#### 4. To help with cleaning, this unit includes a built-in timer.


- 5. We recommend the users to clean their teeth for 2 minutes per cycle. The unit has a 2-minute, auto-shutoff timer.

### Shine Even Brighter...A Few Tips For Best Results!

Clear gel or low abrasion formulations are strongly recommended to enhance the life and performance of your brush head, but pastes may also be used. Abrasive toothpastes, or any toothpaste with baking soda, are not recommended for use with this product. To avoid splashing, switch the power on after placing the brush head in your mouth.

#### 6. Turn on and go slowly, using only light pressure.

Place the bristles perpendicular to flat surface of tooth and slowly guide the brush from tooth to tooth, pausing as you go, so bristles can work into crevices. The OSCILL8™ Rechargeable Toothbrush will automatically provide the correct cleaning action, so there is no need to scrub back and forth. Be sure to run brush along your gum line. Do not apply too much pressure against the teeth. The 8-tuft brush head features tapered bristles oscillating independently to clean effectively without strong pressure. Brush all surfaces, including the outsides, insides, and chewing surfaces. Pay particular attention to back molars.

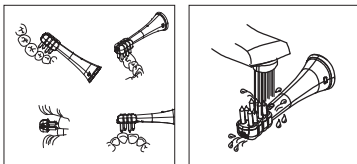
#### 7. IMPORTANT: Rinse the brush head.

After each use, leave the unit on and run warm water through the flushing hole (see picture) and the cleaning ports on front for 15 seconds to thoroughly rinse out toothpaste residue. If not cleaned properly, toothpaste may harden inside the brush head and reduce the life of brush head. Towel-dry power handle and brush head to keep clean and free of buildup.

8. It is highly recommended to flush the brush head per point 7 to prevent toothpaste that could harden inside the brush head and reduce the life of brush head, or cause the eight tufts to stop oscillating or make it hard to start up.

9. After use, turn power handle upside down, shake the handle to remove excess water inside the hollow hole in upper side of handle, drain water out after each use.

10. After use, remove the brush head and any excess water on top of toothbrush body. To prevent excessive water storage in this hollow hole, shake the handle to drain water out after each use. This will prevent water being stored and jamming the waterproof O-ring from moving.
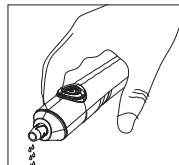
### Limited 1-Year Warranty (U.S. and Canada Only)

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed below, together with your purchase receipt.

California residents need only provide proof of purchase and should call 1-800-334-4031 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

**ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY**

Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

**IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER**

Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

**Service Centers**  
Conair Corporation  
Service Department  
7475 N. Glen Harbor Blvd.  
Glendale, AZ 85307


Please register this product at [www.conair.com/registration](http://www.conair.com/registration)

©2020 Conair Corporation  
Made in China. For more information about OSCILL8™ products, visit [www.OSCILL8.com](http://www.OSCILL8.com) or call 1-800-334-4031

19IN058509  
IB-16083

**OSCILL8™**  
SHINE BRIGHTER

RECHARGEABLE TOOTHBRUSH
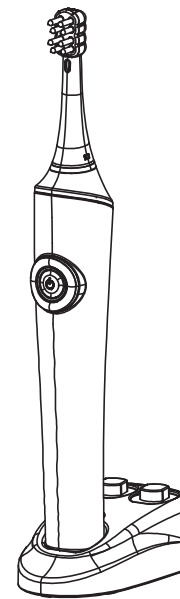
### Instruction Manual • Model NT18

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

Questions or comments? Call 1-800-334-4031  
Visit us at our website: [www.OSCILL8.com](http://www.OSCILL8.com)

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Parte responsable: Conair Corporation, 1 Cummings